



Give your patients more reasons to smile

Opalescence

Tooth Whitening

A brighter smile can give your patients the confidence to smile more. And with Opalescence Go whitening, patients can whiten their teeth when they want, where they want. That's something to smile about.

Powerful, Professional Whitening to Go



- Unique tray material offers a remarkably comfortable fit and easily conforms to any patient's smile
- Molar-to-molar coverage ensures the gel comes in contact with more posterior teeth
- UltraFit tray's superior adaptation ensures the maximum amount of gel stays in contact with the teeth during whitening
- 6% hydrogen peroxide (with wear times from 30–60 minutes)

REFRIGERATION REQUIRED

Opalescence Go Patient Kits



Flavor	6%
Mint	4634
Melon	3592

Each kit contains 10 blister packs with 1 tray each for upper and lower 1 x 20 ml (28 g) Opalescence toothpaste

REFRIGERATION REQUIRED

Opalescence Go Patient Kits 6 pk



Flavor	6%
Mint	4639
Melon	3593

of 10 x upper/lower trays per kit
1 x 20 ml (28 g) Opalescence toothpaste

REFRIGERATION REQUIRED

Opalescence Go Mini Kits



Flavor	6%
Mint	4644
Melon	3599

Each kit contains 4 blister packs with 1 tray each for upper and lower

REFRIGERATION REQUIRED

Opalescence Go Mini Kits 12 pk



Flavor	6%
Mint	4649
Melon	3600

of 4 x upper/lower trays per kit

Procedure


Carefully read instructions and warnings delivered with products before using them.



Remove product from packaging.
"U" - Upper whitening tray
"L" - Lower whitening tray



Position upper tray on teeth.

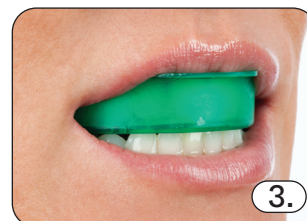

Bite firmly, then suck on tray for 2 seconds.



Remove colored outer tray, leaving the white inner tray on teeth. Repeat the process for the lower tray.



After indicated wear time, remove whitening trays and brush teeth.



Find out more about
Opalescence[™] Whitening
products. Download the
brochure!

Scan QR code

Follow us!



eu.ultradent.blog

© 2025 Ultradent Products, Inc. All rights reserved.

UPP 298-EU 12-2024

ULTRADENT.EU

U2
ULTRADENT
PRODUCTS, INC.